



Complete Solution for OEM Application Testing

The *Transonic Inside™* Evaluation Kit gives you the ability to test the capabilities and performance of the Flowboard and Flowsensor combination that best fits your specific biomedical device/ application. The kit comes complete with the following:

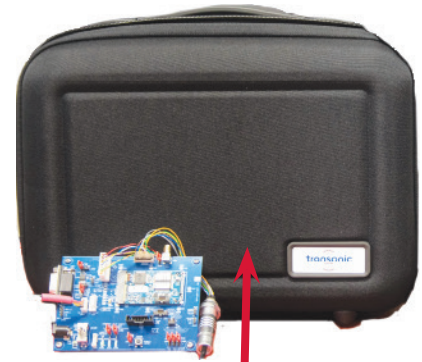

Protective case

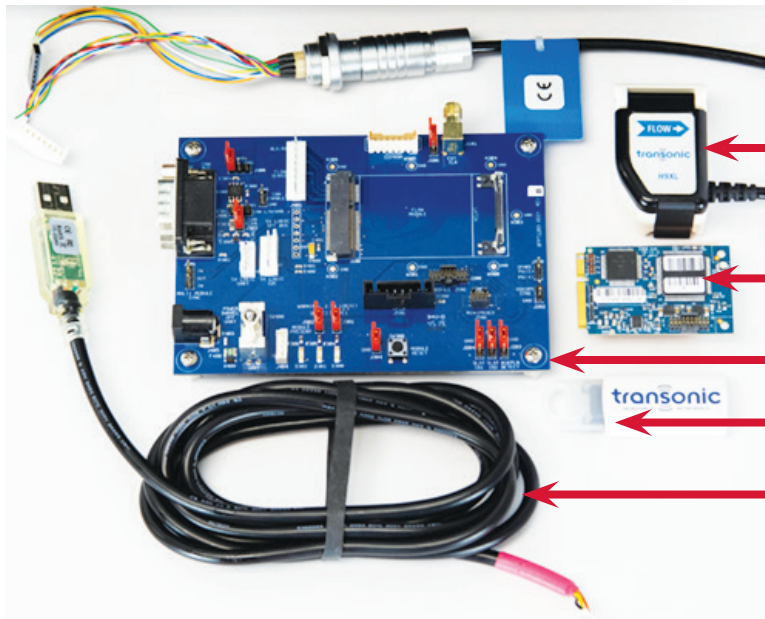

Transonic Flowsensor(s)
with cable assembly

Transonic Flowboard (AFPT1199)

Adapter board


Instructions for use

USB power connector
and adapter

Customizing Your Kit

1. Choose your Transonic Flowboard

Transonic Flowboard		AFPT1168	AFPT1199
Dimensions	cm	4.45 x 1.78 x <.0.64	4.90 x 3.00 x 0.51
	inches	1.75 x 0.70 x <0.25	1.93 x 1.18 x 0.20
Power Requirements		5.0 VDC at 0.25 Watt	5.0 VDC +/- 0.25 at 0.5 Watt
UART Output Interface		SPI, RS-232, UART	SPI, RS- 232, UART, CANopen

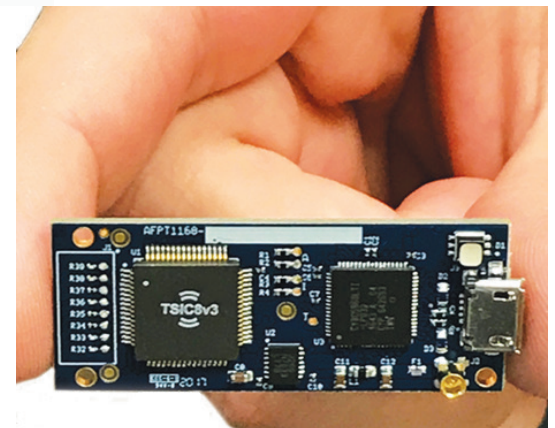

AFPT1168 Flowboard

Bench Testing Solution cont.

Customizing Your Kit cont.

2. Choose your Transonic Flowsensor Size(s)

(Note that while this chart refers to clamp-on sensors, Transonic also provides a full range of inline sensors.)



Tubing Sensor



Open Tubing Sensor

Sensor Size	Tubing ID (Inner Diameter)		Tubing WT (Wall Thickness)		Flow Rate Supported* Liters/min.	
	mm	inches	mm	inches	low flow (1/4 scale)	standard flow scale
2PXL	2.38	3/32	0.79	1/32	0.05 - 0.25	0.1 - 1
3PXL	3.18	1/8	0.79	1/32	0.1 - 0.5	0.2 - 2
4PXL	3.18	1/8	1.59	1/16	0.1 - 0.5	0.2 - 2
5PXL	4.76	3/16	1.59	1/16	0.1 - .0.5	0.2 - 2
6PXL	6.35	1/4	1.59	1/16	0.25 - 1.25	0.5 - 5
7PXL	6.35	1/4	2.38	3/32	0.5 - 2.5	1- 10
8PXL	9.53	3/8	1.59	1/16	0.5 - 2.5	1- 10
	7.94	5/16	2.38	3/32		
9PXL	9.53	3/8	2.38	3/32	0.5 - 2.5	1- 10
10PXL	12.7	1/2	1.59	1/16	1 - 5	2 - 20
11PXL	12.7	1/2	2.38	3/32	1 - 5	2 - 20
12PXL	12.7	1/2	3.18	1/8	1 - 5	2 - 20
14PXL	15.88	5/8	3.18	1/8	2.5 - 12.5	5 - 50
	17.46	11/16	2.38	3/32		
16PXL	19.05	3/4	3.18	1/8	2.5 - 12.5	5 - 50
20PXL	25.40	1	3.18	1/8	5 - 25	10 - 100

* Range includes zero. Any flow peaks exceeding the Max Flow Value (-5 volt to +5 volt range) will be clipped. Measurements below lower linear limit may deviate from the stated accuracy specification. Flow sensors can measure bi-directional flows.

3. Determine Your Calibration Specifications

- Liquid type (e.g., blood, water/saline, Glycerine solution)
- Operational fluid temperature
- Expected flow rate/ range (the chart above includes low flow and standard flow scales).

Since fluid density and temperature alter the transit time of the ultrasound signal and affect the acoustic properties of the tubing, we calibrate the Sensor specifically for these parameters. Operational fluid temperature should be within $\pm 2^{\circ}\text{C}$ of the specified calibration temperature. Specify up to four combinations of tubing, fluid, temperature, and flow rate for calibration. Fluid samples and MSDS information may be required for *custom calibration requests*.

Martine Bosch

Transonic Inside™

Applications Engineer

martine.bosch@transonic.com

John HaberstocK

Transonic Inside™ Market

Development Manager

john.haberstock@transonic.com

transonic
THE MEASURE OF BETTER RESULTS.
info.transonic.com/oem